

# Water Purifier (Distillation)

WS201/221

Production capacity 1.6L/h

Treatment process Distillation


Purified Water Distilled water

## Simple and economical water purifier designed for distilled water production

- Space saving and compact
- By adopting splash prevention design, impurities are prevented from mixing with the distilled water, resulting to a stable water quality
- Designed with empty boiling prevention and overheat prevention functions
- Suitable as washing water for cleaning tools, glassware, etc.



## Treatment Process

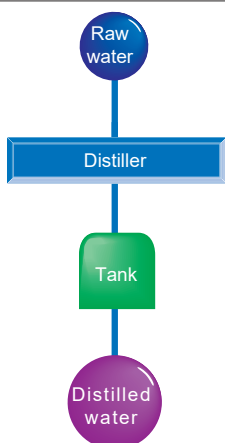

## Specifications

Product code	253061	253062	
Model	WS201	WS221	
Water purifying method	Distillation		
Purified water	Distilled water		
Distilled water production	Approx. 1.6L/h		
Distilled water delivery rate	Approx. 2.2L/min. (When the tank is full)		
Raw water pressure range	1~3×100kPa (1~3kg/cm <sup>2</sup> )		
Safety device	Auto adjustment of cooling water volume		
	Empty boiling / overheat / splash prevention		
Distiller	Boiler	Stainless steel	Hard glass
	Condenser	Hard glass	Hard glass
	Heater	Pipe heater	Built-in quartz glass outer cover
Distilled water storage tank	20L polyethylene tank		
Power source	AC220V 7A		
External dimensions*1	W500×D400×H974mm		
Weight	Approx. 30kg		

● Use water pressure between 1 to 3×100 kPa (1 to 3 kg/cm<sup>2</sup>).

\*1 External dimensions do not include projections.

The length of the power cord is about 2m outside the unit.

## Optional items

No.	Description	Model No.	Product code
①	Stand	AS22	253176
②	Water supply unit	OWH10	253686
③	Drain trap	OWI11	253221



Water Supply Unit (OWH10)



Sintech Scientific (S.E.A.) Pte Ltd

No. 194 Pandan Loop #04-11

Pantech Business Hub

Singapore 128383

Tel : 65-67790007

Fax : 63520877, 63520885

Email : [sales@sintechscientific.com](mailto:sales@sintechscientific.com)

Website : [www.sintechscientific.com](http://www.sintechscientific.com)